



iGo

Portable oxygen concentrator

 **DeVilbiss**
HEALTHCARE



DeVilbiss Healthcare introduces the iGo portable oxygen concentrator. Small and lightweight, yet capable of delivering up to 3lpm of constant flow oxygen meeting the oxygen requirements for most users.

An ideal choice for travelling, the iGo can be operated from a standard 12 V car charger socket or from its lithium ion battery which will give the user up to 5.4 hours of usage (when set at 1 lpm PulseDose mode). The iGo comes complete with a battery, AC power adapter, mains lead and DC power adapter. The iGo can be transported using a wheeled carry case for added portability. The iGo includes the DeVilbiss Oxygen Sensing Device and service alerts to give piece of mind that the iGo is delivering exactly what is required.

A simple, intuitive control panel makes the iGo easy to use and understand. The iGo is robustly designed to deal with the rigours of today's travel. Other than cleaning the gross particle filter, preventative maintenance is reduced as the iGo does not require manufacturer's servicing for the first three years.

Main Features

- Up to 3 lpm constant flow
- Up to 6 lpm pulsed flow (equivalent delivery)
- Long battery life
- Oxygen Sensing Device and audible and visual alarms
- Robust design
- Wheeled carry case available
- 12 V car adapter

Typical Battery Operation Times


Flow Rate	Cont. Flow	PulseDose (20 BPM)
1.0	4.0 hrs	5.4 hrs
2.0	2.4 hrs	4.7 hrs
3.0	1.6 hrs	4.0 hrs
4.0	—	3.5 hrs
5.0	—	3.2 hrs
6.0	—	3.0 hrs



iGo and wheeled carry case



Wheeled carry case showing extendable handle

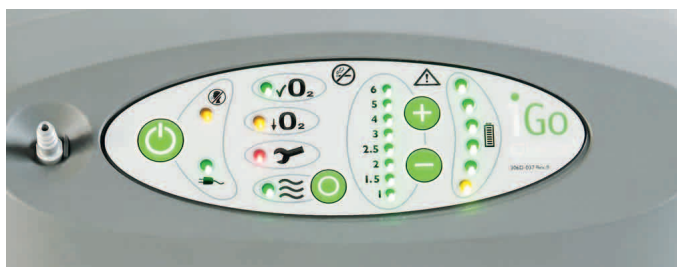

iGo Control Panel



iGo Lithium ion Battery

SPECIFICATIONS

Oxygen Purity	91% ± 3%	Maximum Altitude	0 - 4000 metres
Maximum Constant Flow	3 lpm	Nominal Sound Level	40 dB(A) at 3.0 lpm in PulseDose delivery mode
Maximum PulseDose Flow	6 lpm (equivalent delivery)	Alarm Functions	Low oxygen, No flow, No breath w/failsafe, Low battery, No power, Unit malfunction and High temperature
Calculated Bolus Size	14 cc per l/min		
Weight with/without battery	8.6 kg / 7.0 kg		
Dimensions (H x W x D)	49 x 31.2 x 18 cm		

Distributor:



Sunrise Medical Limited
High Street Wollaston/Stourbridge
West Midlands DY8 4PS
England

Phone + 44 (0) 13 84 44 68 58
Fax + 44 (0) 13 84 44 66 28
www.devilbisshealthcare.com
enquiries@devilbisshc.com