



Discover Life
Again
iGo[®]2



THE GENIUS BEHIND
SmartDose[®]
TECHNOLOGY



Differentiation

Durable Design

Flexible Oxygen Dosing

High Oxygen Output

High Quality Trigger Sensitivity

User Friendly

Durable Design



- Integrated Protective Over-mold
- Recessed O2 Outlet
- Protected Recessed Display
- Uses Pressure Vacuum Technology (PVSA)
- Improved Sieve Bed Efficiency

FLEXIBILITY IN DOSING

SmartDose[®]
TECHNOLOGY




- Patented SmartDose Technology
- Automatically Adjusts to Patient's Oxygen Demand
- Monitors Breathing Patterns and Respiratory Rate
- Dual Setting: Pulse Dose and SmartDose
- Adapts Oxygen Dose to Patient Activity

HIGH OXYGEN OUTPUT



SETTING	Fixed Minute Volume Per Setting	Bolus Volume @ 10 bpm	Bolus Volume @ 20 bpm	Bolus Volume @ 25 bpm
1	260 ml	26,0	13.0 ml	10,4 ml
2	440 ml	44,0	22.0 ml	17,6 ml
3	725 ml	72,5	36.3 ml	29,0 ml
4	880 ml	88,0	44.0 ml	35,2 ml
5	1014 ml	101,4	50.7 ml	40,6 ml

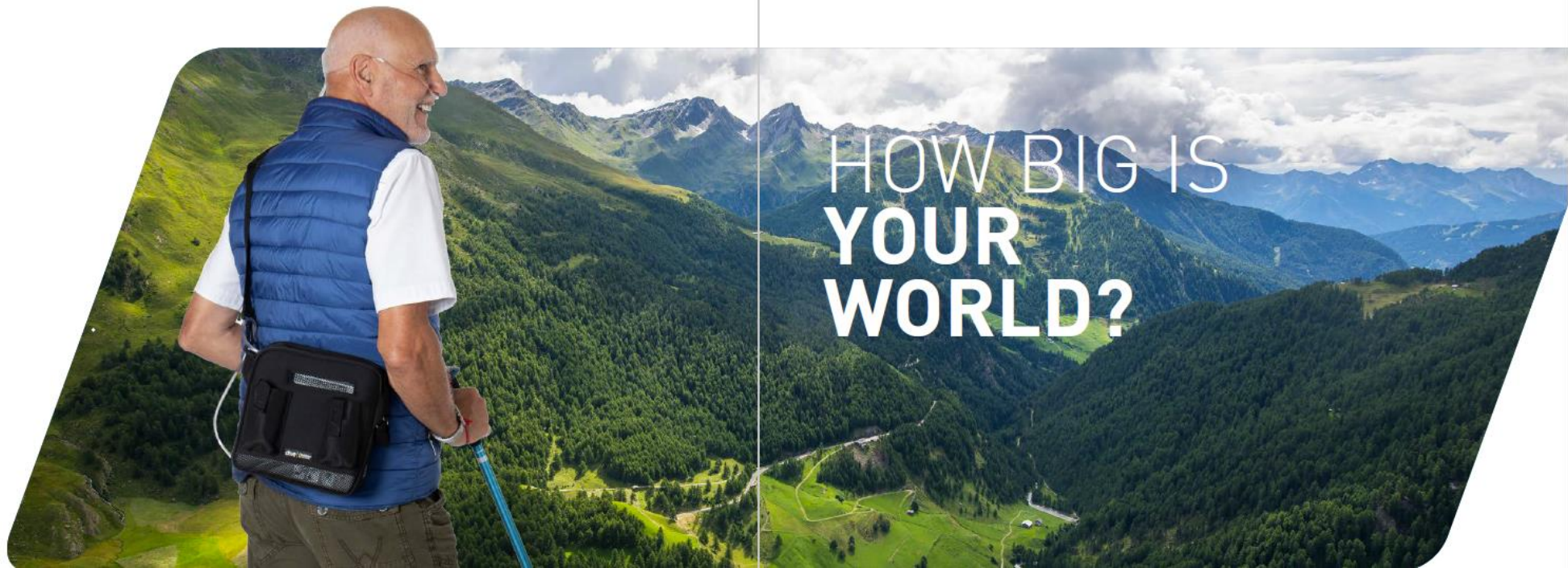

HOW BIG IS
**YOUR
WORLD?**

Advanced
Trigger
Sensitivity
0.05 cmH₂O



DISCOVER
**LIFE
AGAIN**

Ease of Use

Large LCD Display.
Easily accessible even when
used in the carry case.

Easy battery removal.
The iGo2 battery pack can be
easily exchanged without
removing the unit from the
carry case.



Simple to Use Interface

- Vertical Patient Facing Screen
- Protected Recessed LCD Display
- Battery Status Indicator
- Alarm Indicators
- +/- Setting Change
- ON/OFF Button



Accessories



Battery Charger



Case



AC Cord



DC Cord



Carry Case Strap



Battery

Warranty



Key Specifications



Weight	2,2 Kg
Dimensions	21,3 cm x 8,9 cm x 21,8 cm
Max Output	1014 ml/min
O2 Concentration	90% +4/-3
O2 Delivery Settings	1-5 Pulse, 1d-4d/5 SmartDose
Max Breath Rate	40 bpm
Sound Level	< 37 dBa
Altitude	3.048 m above sea level
Operating Temperature	+5 °C to +35 °C
Humidity	15 to 93%
Battery Duration	Min. 3.5 hours @ setting of 2

Key Takeaways

Durable Design

Flexible Oxygen Dosing

High Oxygen Output

Advanced Trigger Sensitivity

User Friendly





drive |  **DeVilbiss**
HEALTHCARE

www.igo2poc.com